



Figure 2

VecTest™

MALARIA Sporozoite Antigen Panel Assay

- *Plasmodium falciparum*
- *Plasmodium vivax* 210
- *Plasmodium vivax* 247

Store at room temperature
FOR INVESTIGATIONAL USE ONLY

NAVIX

Medical Analysis Systems, Inc.
524 Flynn Road
Camarillo, CA 93012

KIT COMPONENTS:

INTERPRETATION OF RESULTS:
POSITIVE FOR DIFFERENT ANALYTES:

Analyte	Control Line	Pf Line	Pv 210 Line	Pv 247 Line
Negative	+	-	-	-
Pf	+	+	-	-
Pv 210	+	-	+	-
Pv 247	+	-	-	+

• The test should be interpreted positive only if at least two lines appear on the strip.

• The test is negative when only one line (control line) develops on the strip.

• Discard strip if no line develops. The test results are not valid if control line does not develop.

PROCEDURE:

PREPARE MOSQUITOES FOR TEST → HOMOGENIZE SAMPLE → REMOVE THE GRINDER, PLACE THE DIPSTICK IN THE SAMPLE TUBE AS SHOWN IN THE FIGURE

PLACE MOSQUITO PARTS (HEAD AND THORAX) INTO THE GIVEN TUBES WITH 250µL (13 DROPS) OF GRINDING SOLUTION

WAIT 15 MINUTES BEFORE READING RESULTS

Signal pattern:


Control —
Pf —
Pv 210 —
Pv 247 —

Place your test strip in the box for comparison



REPLACEMENT SHEET


Figure 3

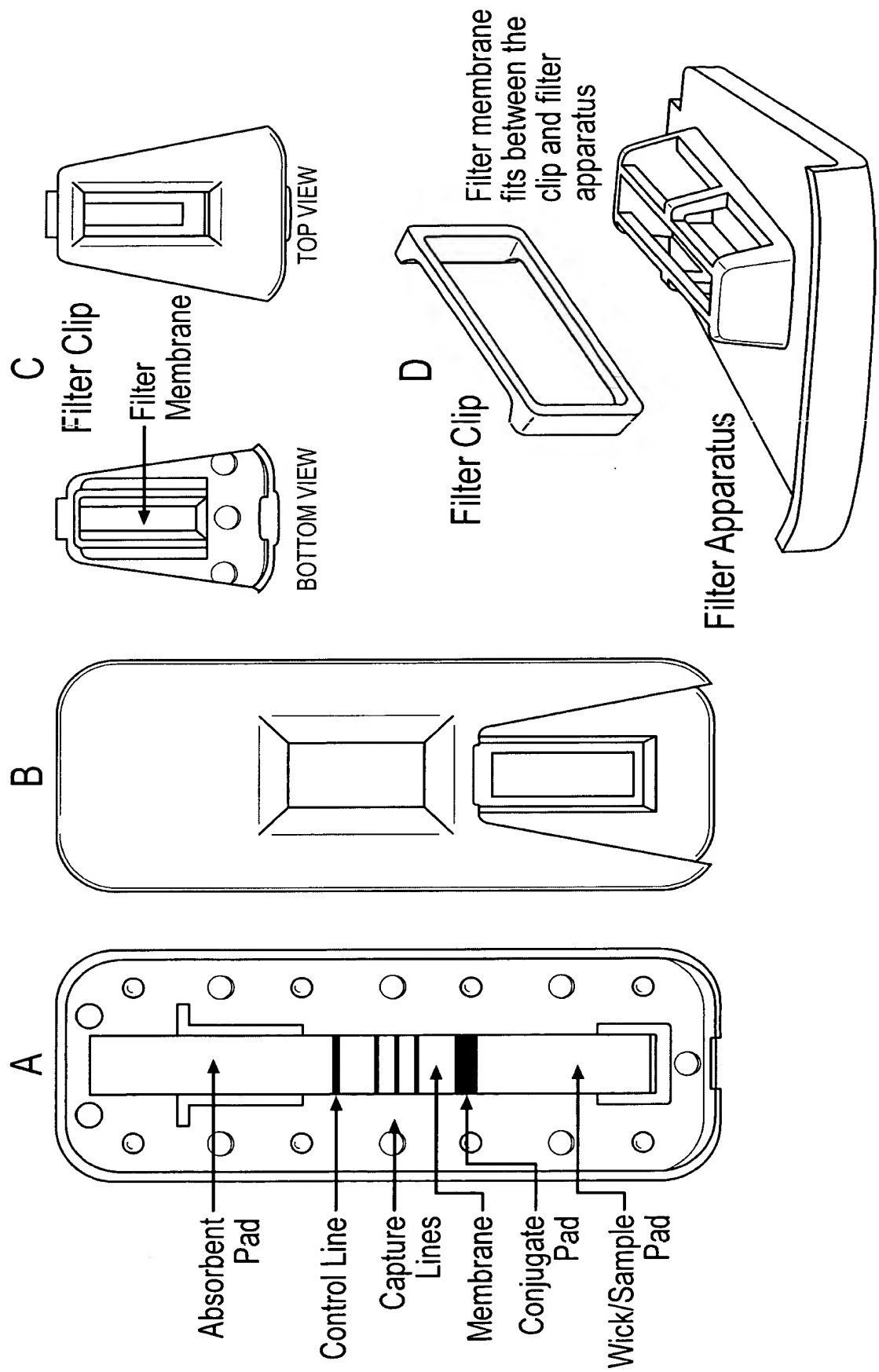



REPLACEMENT SHEET


Plasmodium falciparum CS Protein Detection
ELISA versus Dipstick: Figure 4

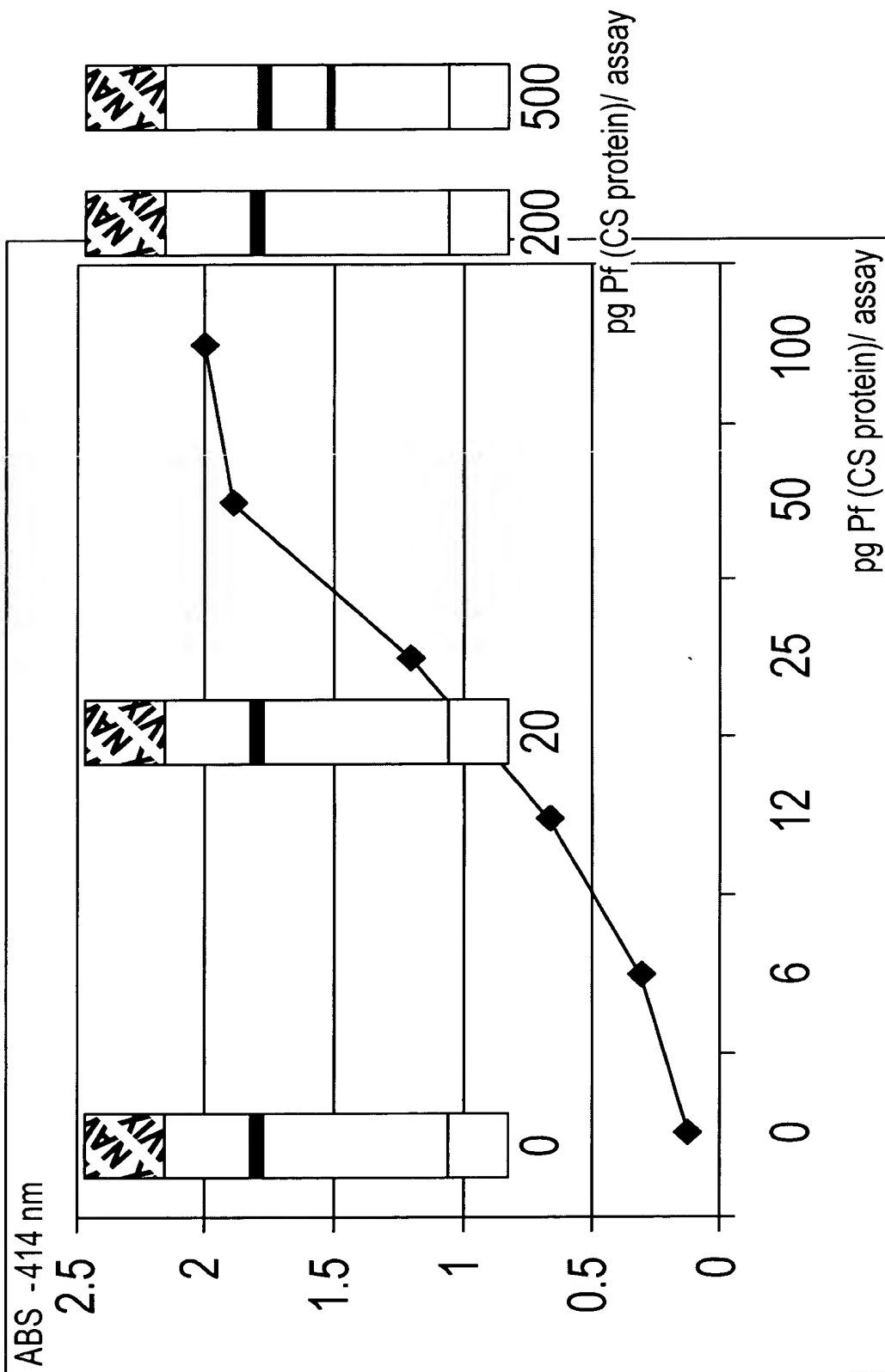




REPLACEMENT SHEET

Plasmodium vivax 210 CS Protein Detection
ELISA versus Dipstick: Figure 5

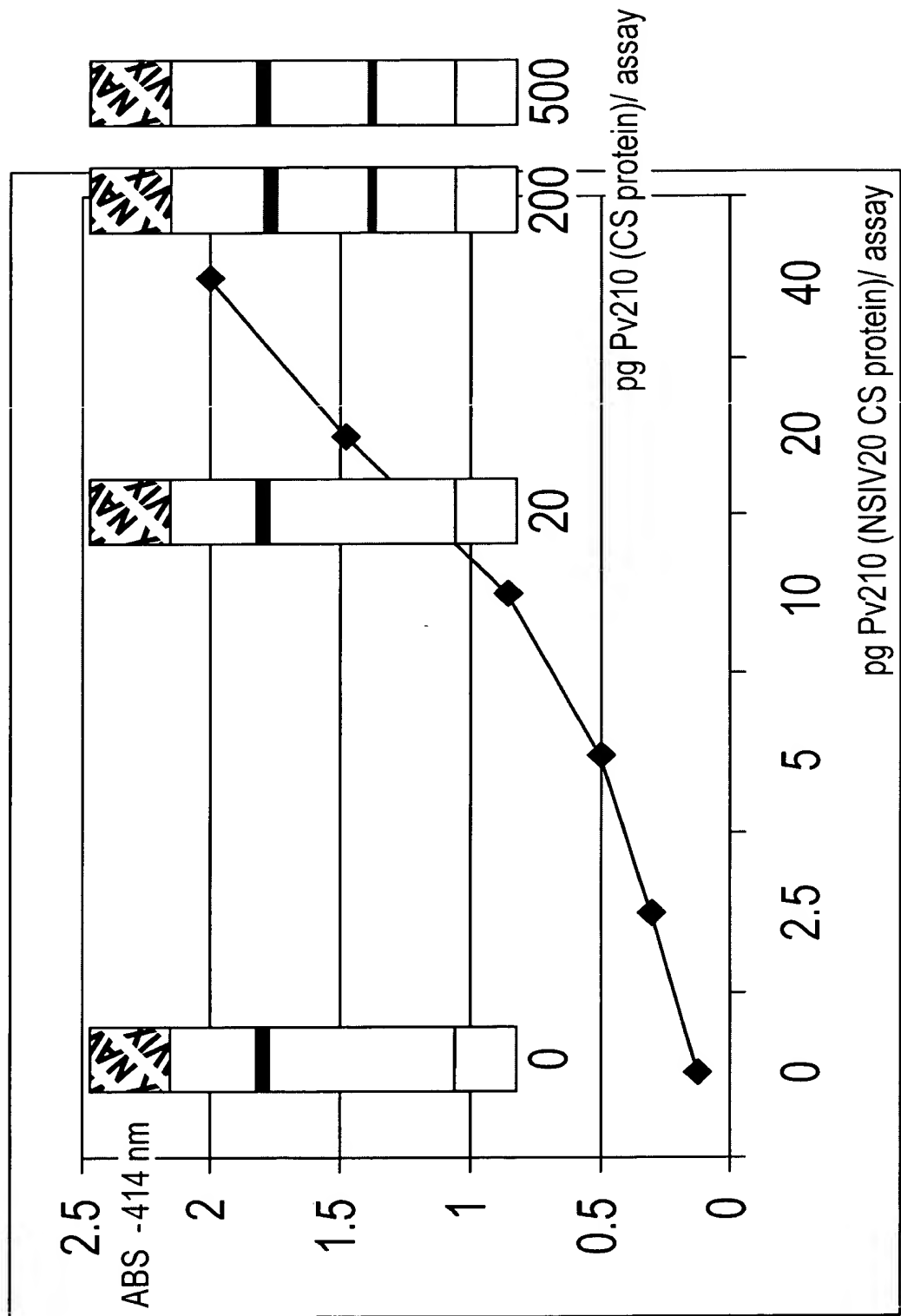

Plasmodium vivax 247 CS Protein Detection
 ELISA versus Dipstick: Figure 6

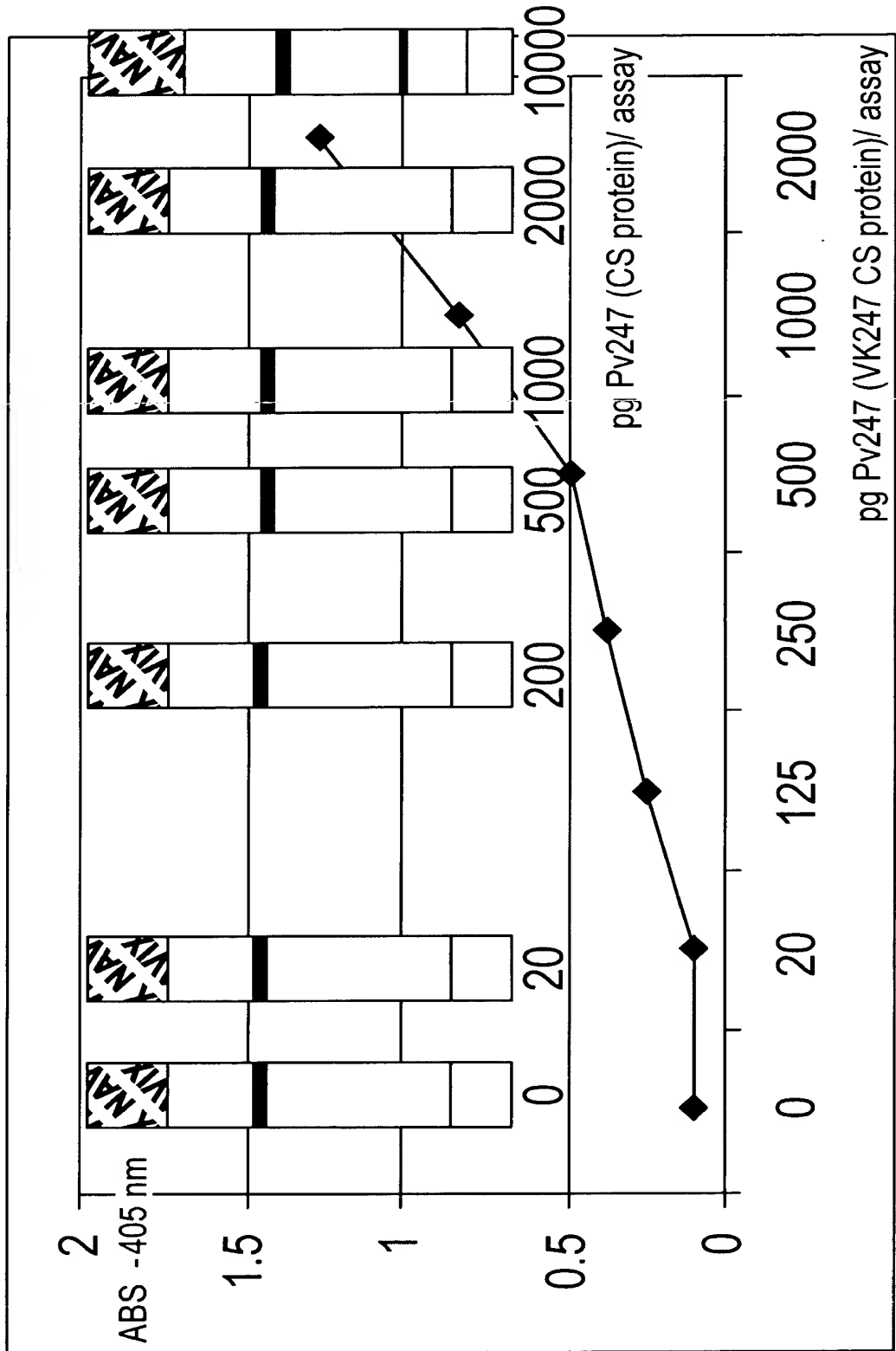



Panel Dipstick Detection of *Plasmodium* Sporozoites in Infected Mosquitoes: Figure 7

Mosquito (MQ) extract:

Pf infected mosquitoes Pv210 infected mosquitoes Pv247 infected mosquitoes

1MQ, 1MQ, 2MQ, 3MQ, 5MQ

1MQ, 1MQ, 2MQ, 3MQ, 5MQ

1MQ, 3MQ, 5MQ,



>400

300

>400

30

300

0

0

400

35

75

75

12.5

<12.5

50

Pf sporozoites
(ELISA result)

Pv210 sporozoites
(ELISA result)

Pv247 sporozoites
(ELISA result)

Sensitivity in ELSA:

~25-50 Pf sporozoites

~25-50 Pv210 sporozoites

~50 Pv247 sporozoites

REPLACEMENT SHEET